



**Mirasoft Retractors**

These soft & pliable retractors feature long handles allowing maximum patient comfort. Autoclavable.

Adult & Child Sizes- Set of 2



**Soft-Dex Retractors**

Provides a full view of the anterior & posterior regions. Autoclavable.

Adult & Child Sizes- Set of 2



**Mirahold Retractors**

Clear plastic with long handles. Can be cold sterilized. Autoclavable.

Adult & Child Sizes- Set of 2



**Spandex Cheek Retractors**

The original retractors. Adult, Child & Teen Sizes- Set of 2



**Retract-Ezz**

- Allows unrestricted access.
- Great for occlusal photos & quick direct bond repair.
- Extremely comfortable for patients.
- Angle head helps retract lips.
- Cold or Heat autoclave up to 290°F.



**Fog Free Intraoral Photography Mirror System**


- Micro fan in base prevents fogging
- 3 L.E.D. lights provides direct lighting
- Cordless – charger included
- 4" x 2.5" x 1" base serves as a handle during photography



Kit includes: base, #3 mirror, charger & barrier sleeves



**Available Mirrors for Fog Free System**

Interchangeable highly polished stainless steel mirrors with a special coating process

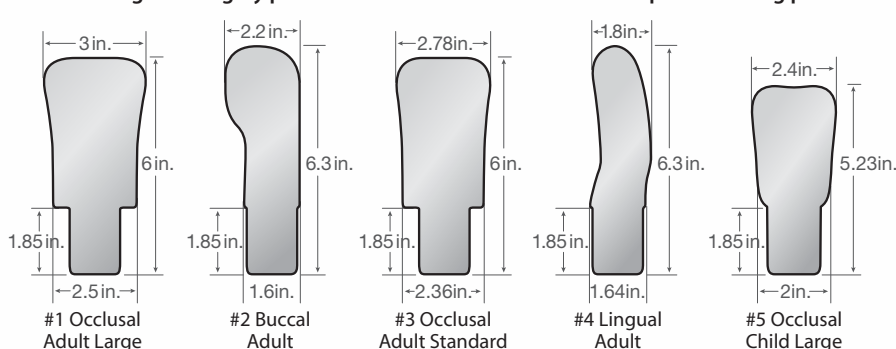


**Glass Photo Mirrors**

- Autoclavable
- Double Sided
- Rhodium Coated
- Front Surfaced
- Available in 5 styles
- Made with Lapidated Crystal

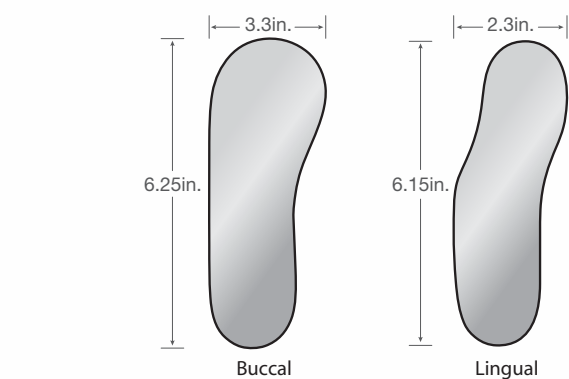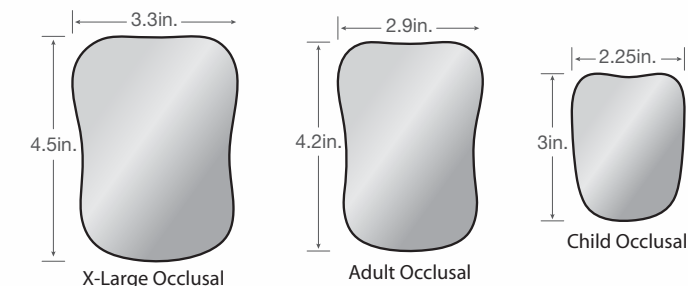

**Stainless Steel Photo Mirrors**

- Highly Reflective Surface
- Precision Laser Cut
- Available in 6 styles

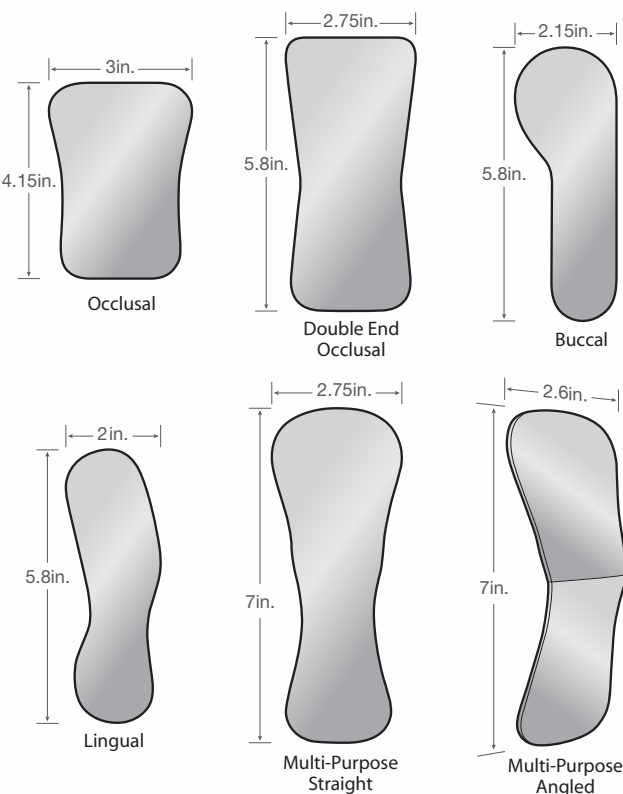

**Polar Blast™ Freeze Spray**

Helps to assure maximum adaptation of trays, retainers, mouthguard materials and other plastics by stopping all thermal action.

Increase flexibility on thermally activated arch wires. Works great with heat activated wires by making wires more pliable.

10 oz. cans contain 4 oz. more than most other name brand freeze sprays.



**Endo Freeze™ Refrigerant Spray**

Spearmint Scented Endo Freeze is a quick acting spray specially formulated for testing pulpal vitality. The low liquid temperature of -15.7°F (-26.2°C) insures immediate results from vital teeth.

Safe for use in any environment, Endo Freeze™ is non-flammable and non-explosive. Ideal for on-the-spot pulp testing.

Select Dental Endo Freeze™ contains zero CFC and is environmentally safe for the Ozone.

6 oz. cans.



**Contouring Kit**

Removes interproximal tooth mass to produce additional intra-arch space. May also be used for opening contact points, recontouring or removing overhanging fillings. Reusable aluminum-oxide contouring blades.

Kit Includes 1 Handle and 12 Assorted Strips

**Additional Strips Available:**

- Single sided strips- package of 10 pcs.
- Double sided strips- package of 10 pcs.

**Strips Available In:**

Fine, Medium and Coarse



**Diamond Disc Guard**

The stainless steel Diamond Disc Guard protects the mouth from spinning Diamond Discs. Fits all 7/8" diamond discs.

For use on doriot style handpieces.



**Enamel Hand Stripper Kit**

Manual operation means more control and increased safety. Can be used up to the first molar and accepts all sizes of steel abrasive strips and many types of linen, cloth or plastic strips. Single and double sided.

1 abrasive strip can be used for several applications. Using the correct width prevents gingival tissue injury.

Kit Includes 1 Stainless Steel Enamel Hand Stripper and Abrasive Strips Sample Pack.

**Additional Strips Available:**

- Single sided strips- package of 10 pcs.
- Double sided strips- package of 10 pcs.

**Strips Available In:**

Fine, Medium and Coarse

**Widths Available:**

4mm, 6mm, 8mm



**Mini Stripper**

- Small and easy to handle
- Complete with 16 refill strips (single side only)

**Additional Strips Available:**

- Serrated strips- package of 10 pcs.
- Steel strips- package of 10 pcs.
- Diamond strips- package of 10 pcs



**Cello Stripper**

The Cello Stripper Holder accepts 2 sizes of steel abrasive strips- 2mm and 3mm; single-sided and double-sided.

Kit Contains:

- Stainless Steel Strip Holder
- Sample Package of 10 Abrasive Strips (2mm or 3mm width, single-sided or double-sided with and without blanks)

**Select Dental Model Forming Bases** give you great looking models without the hassle and mess associated with model trimming. These tough yet flexible bases offer easy separation as well as long life. All molds can also be purchased individually.



**Rubber Model Forming Bases**

Ortho Molds- regular set of 6 pieces Available in Upper or Lower (Sm, Med & Lg) Approx. 1/2" high in front and 5/8" high in back



**Universal Molds**

Can be used for Upper or Lower Models Approx. 1/4" high in front and 1/8" high in back



**Mixing Bowls**

Flexible, chemical resistant, neoprene rubber mixing bowls. Perfect for mixing plaster, dental stone, and silicone casting materials.

- Small- 3 1/4" W x 2 1/4" H 150 ml
- Medium- 4 1/4" W x 2 3/4" H 350 ml
- Large- 4 3/4" W x 3 3/8" H 600 ml
- X-Large- 5" W x 4 1/2" H 850 ml







### Megatray™ Light Cure Custom Tray Material

Megatray™ is a light cure tray material supplied in economical prefabricated sheets- eliminating mixing and mess.

Megatray™ features good working time along with excellent malleability and adaptation to model. Easy trimming characteristics reduce error as well as effort.

Megatray™ is odorless and uniform in thickness with low shrinkage and exceptional strength.

2mm- Available in Blue and Pink- 50 Trays per box  
1.4mm- Available in Pink- 50 Trays per box



### Maxi-Cure™ Curing Light

The heavy-duty construction makes the MAXI-CURE™ light cure box a reliable way to cure materials.

Light uses four 9 watt bulbs (included) and cures most light activated composites including Megatray.

Complete with a 5- and 10- minute timer, this unit measures 11" x 5" x 4 1/4" .

**Full One Year Warranty**



### Vibrator Variable Speed

Allows the technician to vibrate a mixing bowl or flask at the correct speed for the plaster, stone or viscous solution to de-bubble.

The rubber top is removable for easy cleaning and replacement.

Four soft rubber suction cup feet prevent spills caused by creep with a firm hold even during the heaviest vibration cycles.



### Personal Ultrasonic Cleaner with Detachable Tank Model 3900

- 5 minute Auto Shut Off
- Steel Tank Capacity: 0.3pt (140 ml)
- Unit weight: 1.5 lbs.
- Dimensions: 4 3/4" w x 4 3/4" h x 7" d



### Digital Ultrasonic Cleaner Model 3800(A)

- Steel Tank Capacity: 1 pt (600 ml)
- Unit weight: 2.2 lbs.
- Dimensions: 8" w x 5 3/4" h x 5 1/2" d



### Digital Ultrasonic Cleaner Model 4800

- Steel Tank Capacity: 2.4 pt (1,375 ml)
- Unit weight: 5 lbs.
- Dimensions: 7" w x 6" h x 9" d



### Digital Ultrasonic Cleaner with Built-in Heater Model 4820

- Steel Tank Capacity: 5.5 pt (2600 ml)
- Unit weight: 6 lbs.
- Dimensions: 11" w x 7" h x 8 3/4" d



### Articulators

These solid brass or die cast and chrome plated articulators are constructed with stainless steel and brass fittings- sturdy and rugged for years of use, case after high quality case.

The two piece articulator design provides standard lateral and protrusive movements.

The upper and lower parts are easily separated.



### Crucibles

Created using a special manufacturing process & made from the highest quality fine grain quartz, these crucibles are ideal for all types of metal. Even heat distribution provides for finer casting.

- Non-contaminating
  - Long Lasting
  - Withstands high heat
- \*Packaged Individually



### Brass Denture Flask Upper & Lower

Features forged bronze alloy construction.

### 2 Case Compress

For use in curing two Flasks simultaneously.

### Reline Jig

Chrome plated, 2 post jig for accurate reline applications.



### 207 Micromotor

- 35,000 rpm
- Features Forward – Reverse.
- Non Stage Speed System (0 – 35,000rpm)
- Includes variable speed foot control



### Portable Micromotor

- Rugged and rechargeable, yet compact
  - Stepless speed adjustment from 0–25,000rpm with great torque throughout
  - 6 hour charge time with an incredible 8 hour continuous run time even under heavy load
  - Sealed Ni-MH battery assures extremely long life
  - E type handpiece comes with a doriot straight handpiece nose cone and will accept any kind of E type front end
  - Comes complete with holster and charger
  - Other front ends available
- Full One Year Warranty**



### Easy Vac™ Vacuum Former

The Easy-Vac™ Vacuum Former features a built-in vacuum and accommodates square or round material. Sturdy construction ensures years of reliable use.

- Hex Post allows smooth up & down movement
- Reflective heating element casing disperses the heat evenly over any thermoforming material
- Dial Timer auto shut-off prevents overheating
- Compact size with "no-worry" operation makes it easier to store and use in your lab
- Noise Filter provides quieter operation
- 2-Fuse Operation keeps the machine from overloading the circuits

Power Supply- 110 – 220 V AC / 50 – 60 Hz  
Heat Output- 550 W • Motor Output- 1250 W



### FX 45

#### L.E.D. Cordless Curing Light

- Features FAST CURE cycle as well as RAMP, CONSTANT & PULSE CURING cycles
- Built-in Radiometer
- Digital display of curing mode & timer
- Long lasting & replaceable Li-ion battery
- Up to 120 min. run time per charge
- Sterilizable Curing Tip protects L.E.D. Module
- Battery charging station with light intensity meter



### MIGMA 200 Alginate Mixer

#### Your First Step to a Lasting Impression

- Eliminates hand mixing while guaranteeing a homogeneous bubble free impression, using any alginate you choose
- Precision engineered in Switzerland
- Measures a compact 8" x 12" x 10"
- Includes 2 large lidded mixing bowls and one mixing spatula



### Prestige Classic Autoclave

- Light, compact & portable, weighs only 9.9 lbs.
- Fully automatic- single, one step operation
- Large capacity chamber 8" x 9"
- Safe and reliable
- Temperature: 250°F (121°C)
- Sterilization Time: 16 minutes
- Cycle Time: Approx. 30 minutes
- Includes: Cassette Rack and 3 Cassettes



### MR Dual Pin & Sleeve

The metal sleeve has a high luster plating that is abrasion resistant, allowing the pin to be removed and replaced frequently without any change to the relationship of the overall model.

- Extremely well made, accurate & economical
- Virtually eliminates rotation
- Compatible with virtually all pin setting systems
- Precise fitting metal pin and metal sleeve with improved retention
- Die Guard for clean professional models
- Can replace: PX114, PX115, PX153, P-48, Pin Mates & Denerica 346 & 326
- Available in boxes of 100 pcs. & 1000 pcs.



### MR Flats Pin & Sleeve

- New concept in Dowel Pins
- Flat side prevents rotation associated with round pins
- Tapered side assures pins are always replaced properly
- Great for anterior and posterior Dies
- Plastic sleeve has die guard built in (No need to assemble- just snap off end)



### Diamond Saw & Mandrel

- 2 3/4" diamond coated blade with special safety mandrel
  - Cut neat clean dies from casts
  - Use on standard dental lathe
- Set contains: Diamond Saw, Mandrel and Wrench  
(All items also available separately)



### Slim-Line™ View Box

#### Ultra Flat- Only 3/4" Thick

The attractive Slim Line™ View Box is only 3/4" thick and can easily hang on a wall, lie flat or at an angle with built-in leg support.

Perfect for patient presentations, cephalometric tracing as well as day in and day out use.

Generous 12" x 8" viewing area accommodates a large variety of film sizes- the outer dimensions measure 14 1/2" x 10 3/4".

Cold cathoid technology assures reliability as well as extremely long lasting bulb life.

**Unit includes 1 AC Adapter**

**Full One Year Warranty**



### Digital Wax Pot

- Precise digital display of temperature at a glance
- Easy to use "Up" & "Down" temp. control buttons
- Temperature Range: 122°F – 239°F
- Efficient: quickly reaches & maintains desired temperature
- Accurate: makes perfect copings every time
- Compartment Depth: 5/8"
- Size: Approx. 5 1/8" W x 6 1/2" D x 2 5/8" H
- Weight: 2.5lbs
- 115V 60Hz



### Analog Wax Pot

- Temperature dial provides quick control with "Low" and "High" settings
- Temperature Range: 122°F – 239°F
- Efficient: quickly reaches and maintains desired temperature
- Accurate: makes perfect copings every time
- Compartment Depth: 7/8"
- Size: Approx. 5" W x 4 1/16" D x 2 1/4" H
- Weight: 2lbs
- 110V 60Hz



### Laser Pin Setter

Compact and rugged, The Laser Pin setter is compatible with virtually all pin setting drills and dowel pins, from PX 114 & 115, Medium Dowel Pins and of course MR Pins. Its high torque motor makes drilling easy, while its focused laser beam ensures accuracy.



### MR Pindex Style Drills

Fits virtually all pin setting machines. Drill the perfect hole for the MR Pins, #2 Medium Dowel Pins, and Pindex Pins.

Comes in packs of 3 Drills.



### Delivery Bags

- Poly-lined to withstand humidity and wetness.
- Extra thick 60lb white paper for a clean and professional look
- Large gusset for a roomier articulator fit
- Equipped with a tin tie closure for a secure seal
- Available in 2 sizes:  
Regular- (5.5" x 3.5" x 11") Case of 100 or 500 bags  
X-Large- (6" x 4" x 13") Case of 100, 250 or 500 bags



### Brass Dowel Pins

- Available in 3 sizes (Sm, Med & Lg)
- Packaged in units of 100 pieces and 1000 pieces



### Brass Stick Pins

Non-corrosive, high quality brass stick pins ideal for all impression material.  
Packaged in sets of 100 pieces and 500 pieces



### MR Auxiliary Pins

- Nickel plated dowel pins
- Packaged in sets of 1000 pieces



### PinPro

#### Dowel Pin Adhesive

The only Cyanoacrylate Adhesive formulated specifically for dowel pin placement. Bonds virtually everything with its special gap filling formulation.

Works on oily surfaces too!  
1 Fluid Ounce